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USER MANUAL

ENGLISH

06

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21


FRANCAIS

35

PACKAGE DETAILS

**ON YOUR PACKAGE,
YOU WILL FIND :**

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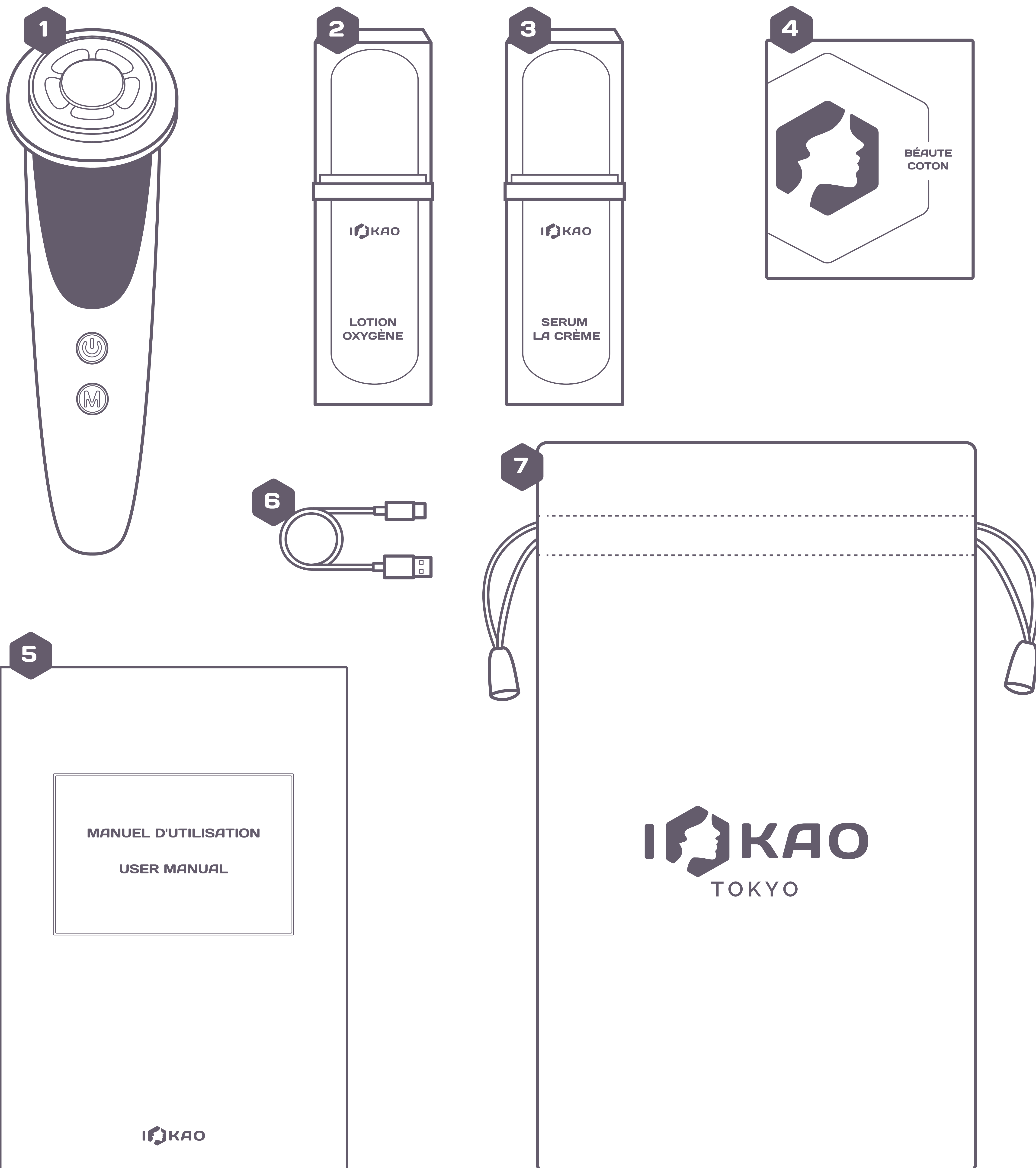
M button

PRECAUTION

Avoid charging the device with external cables.

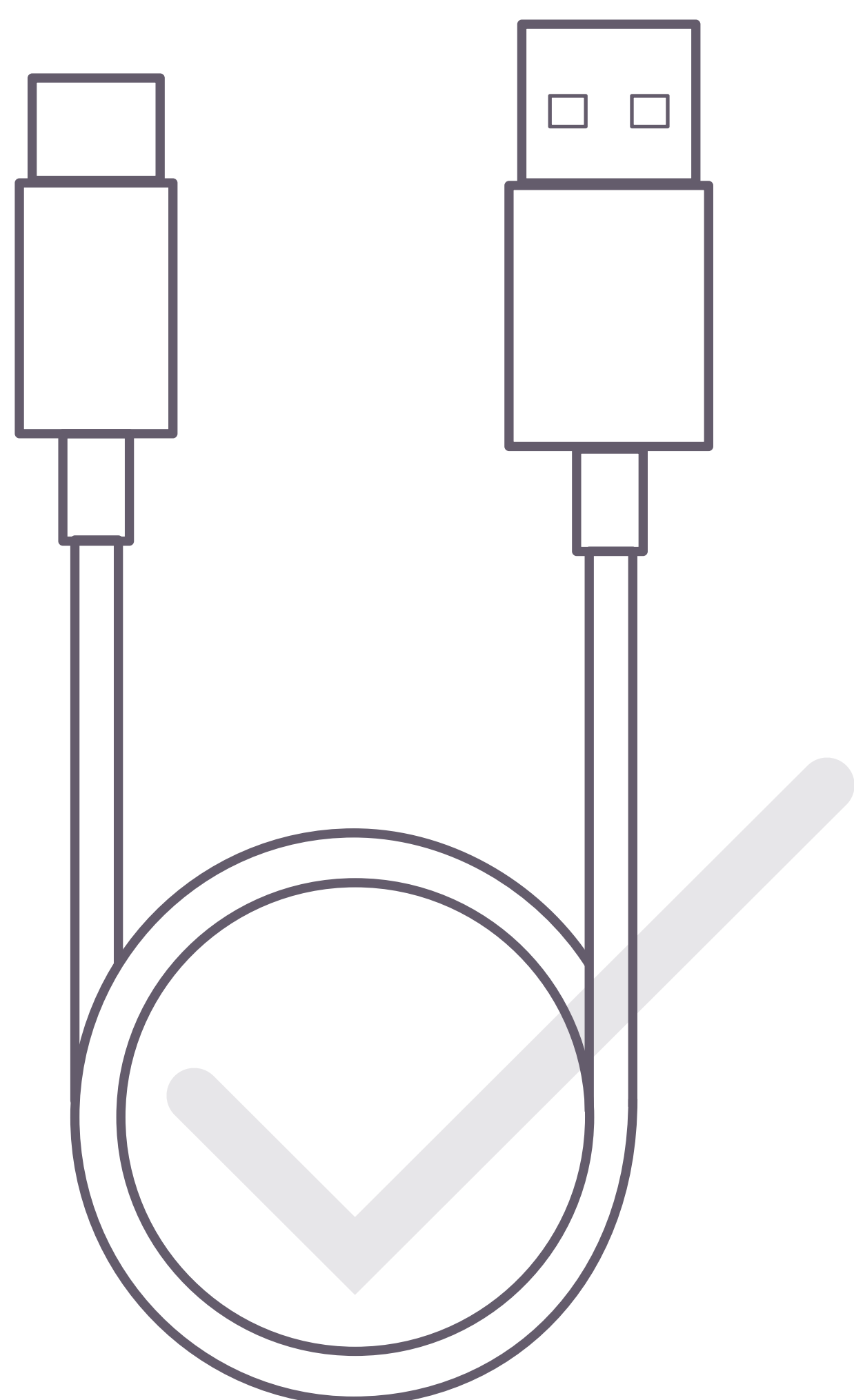
Do not use water to clean the equipment.

OPERATING MODE


1 POWER ON/OFF

Press the on/off button for three seconds to turn on/off.

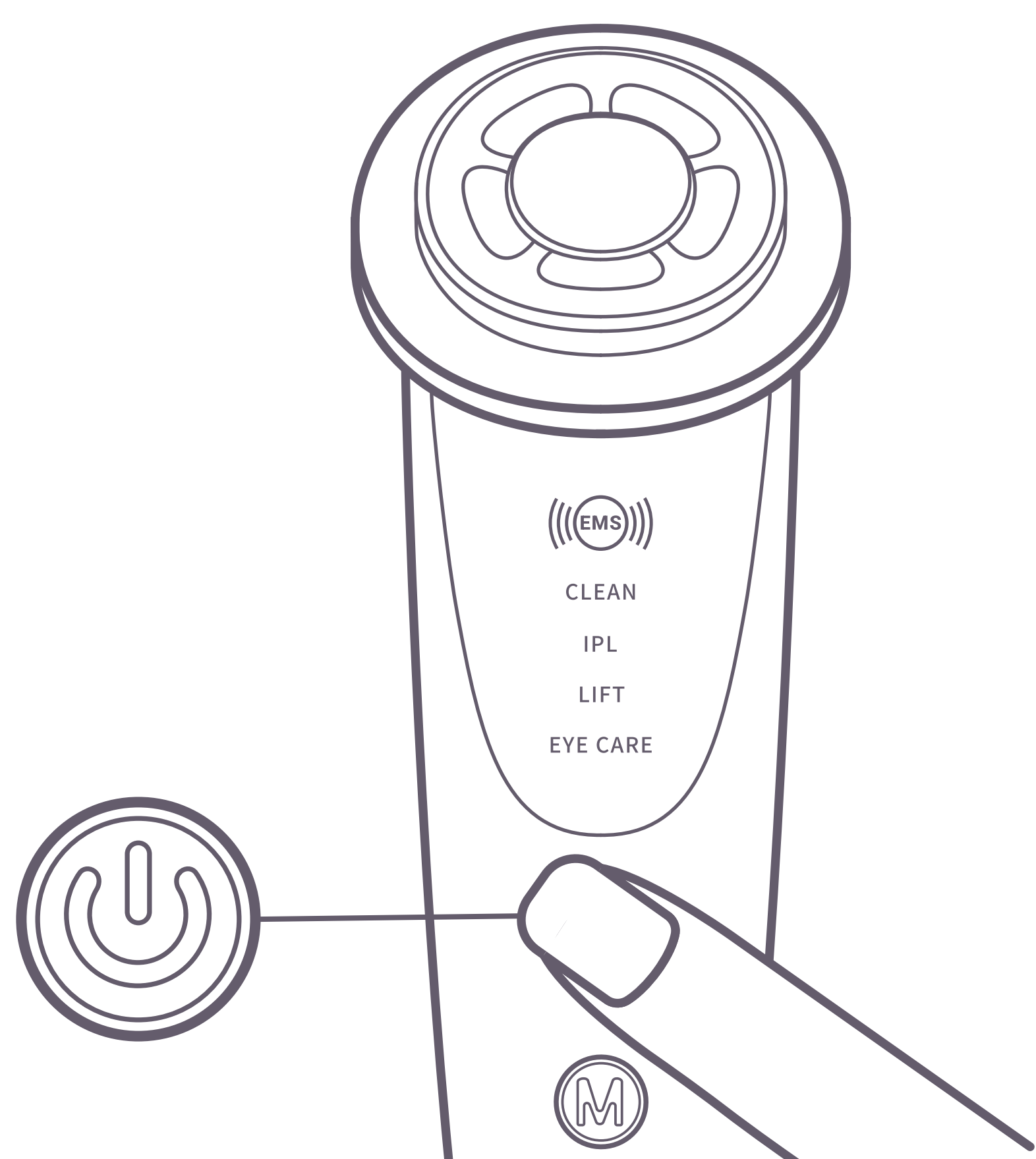
2 SELECT MODE

Briefly hold the on/off button to switch modes (Clean, IPL, Lift, Eye Care)

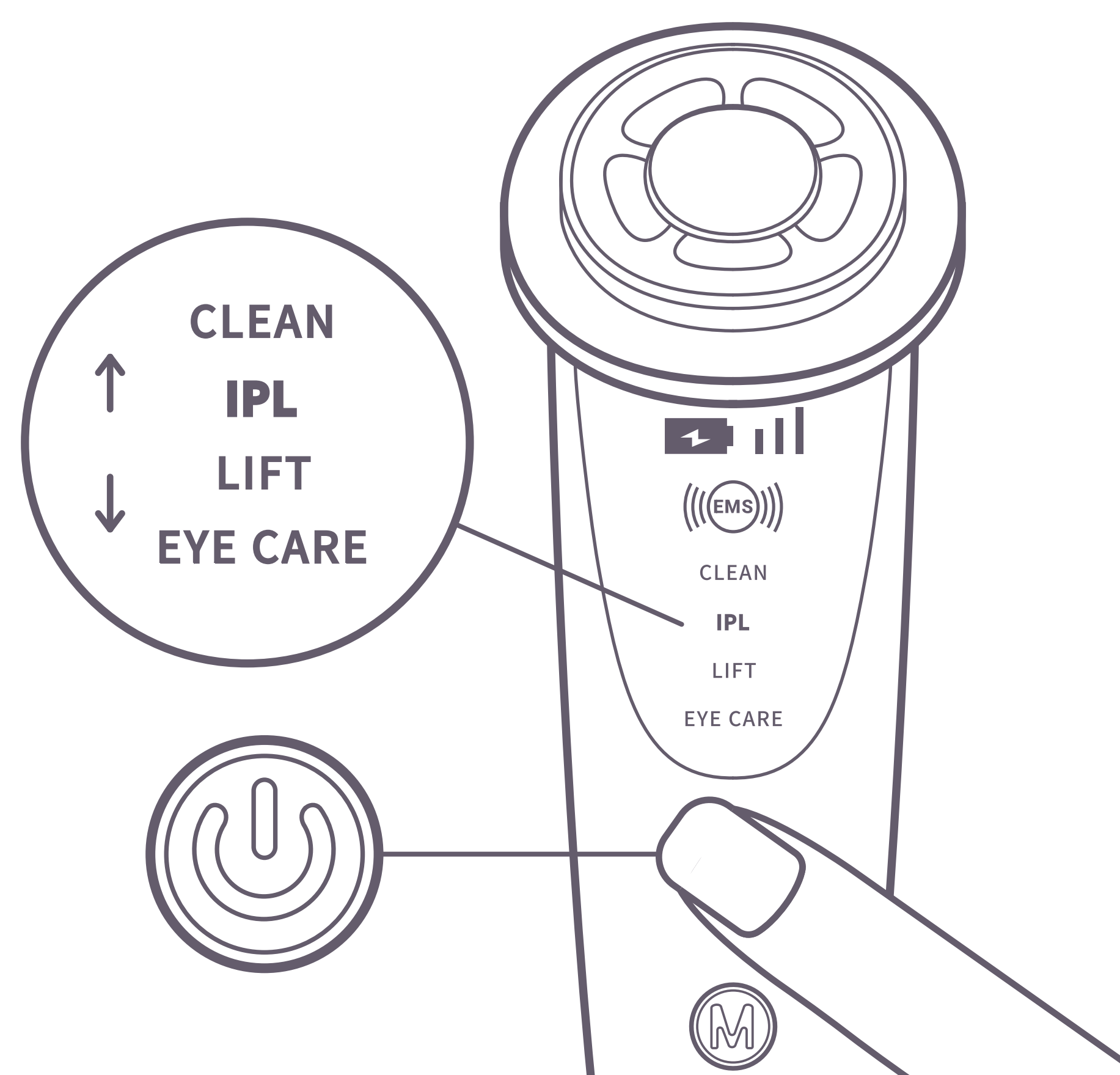
3 SET VIBRATION INTENSITY

By pushing the "M" button, you may select the vibration level (low, medium, or high) for each mode (except Eye Care Mode).

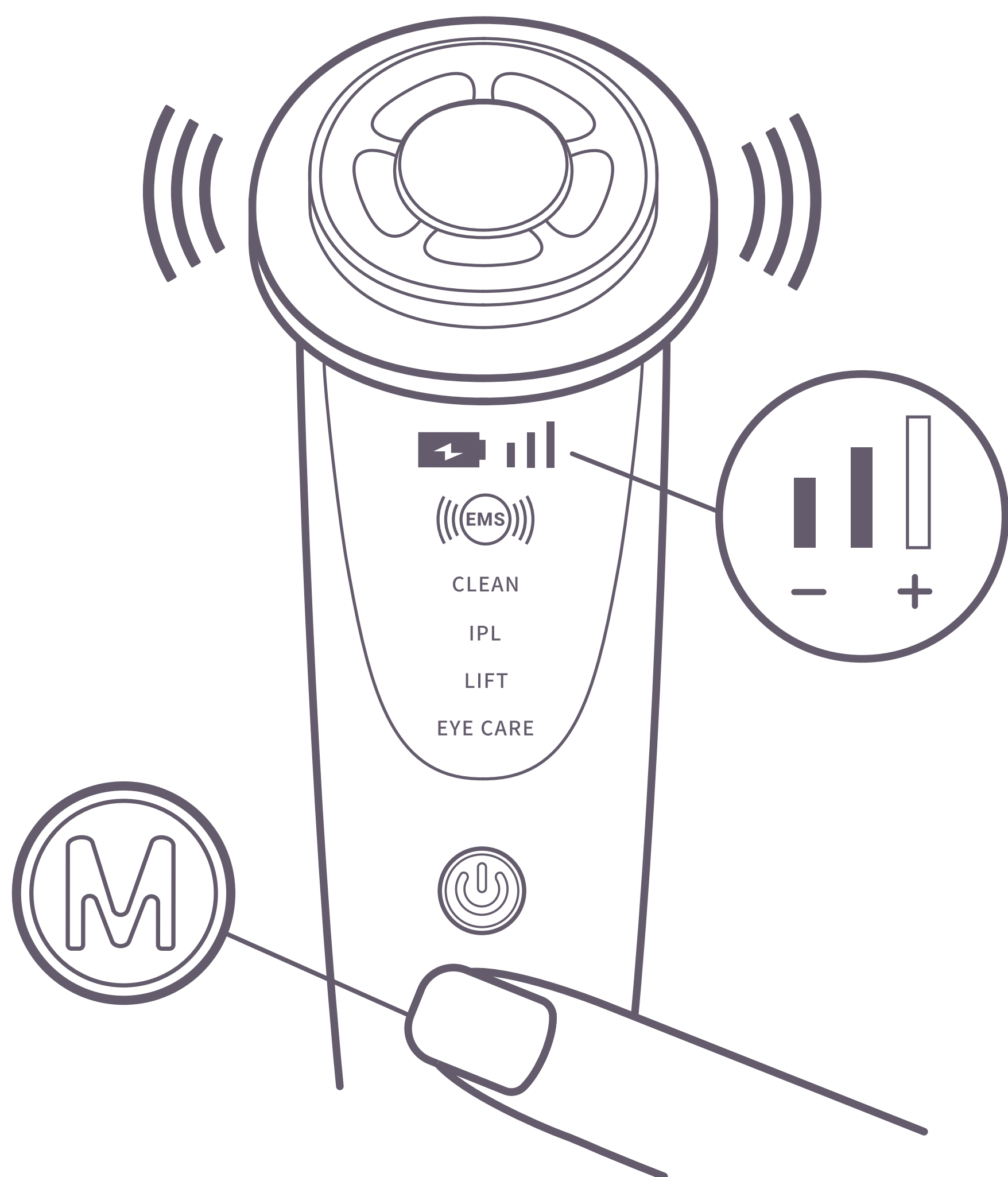4 EYE CARE LIGHT

When you select this mode, you can choose whether to initiate red light therapy or not by pressing the "M" button.

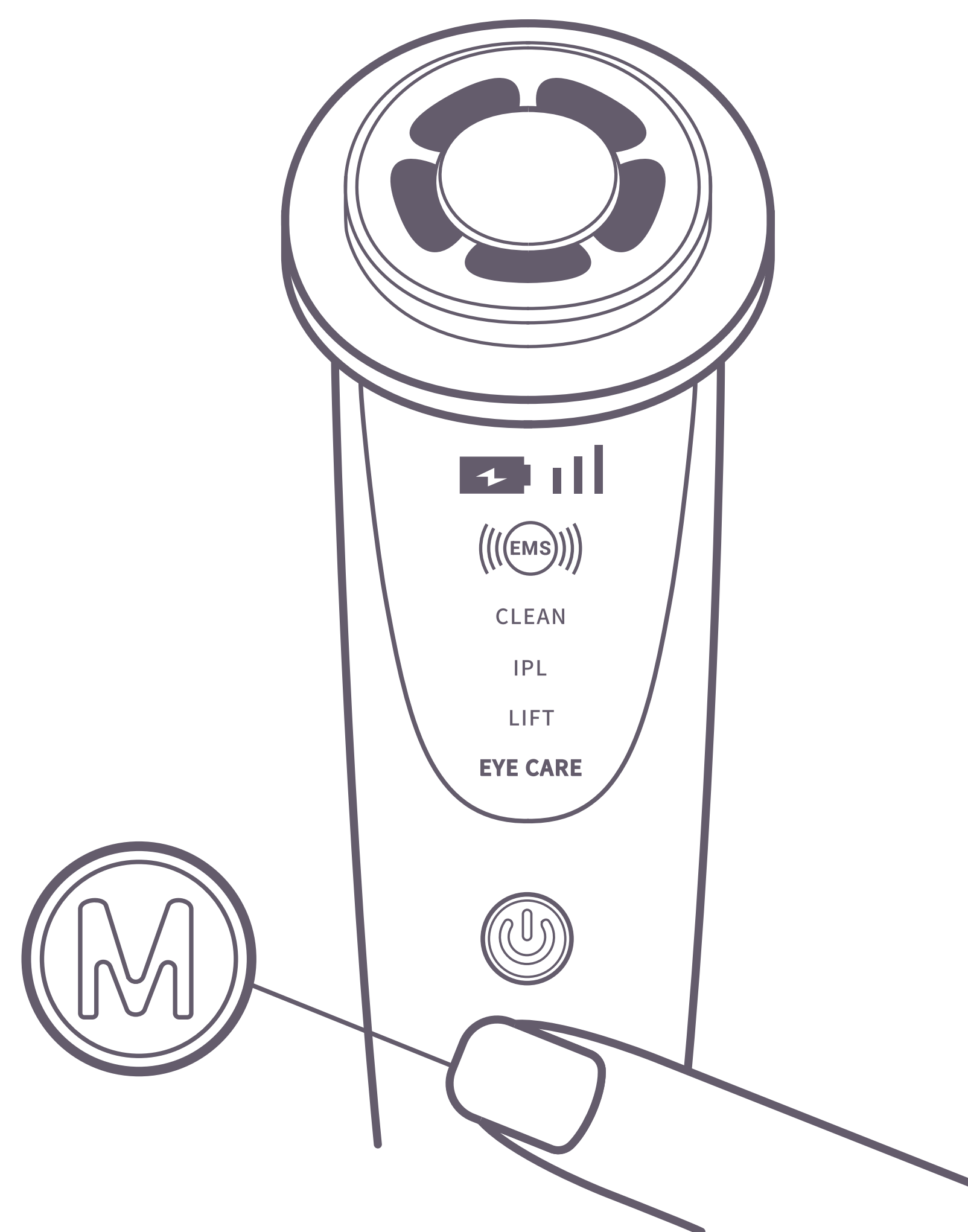MODE CHOICE

CLEAN

Select the Clean mode to remove dirt and stains from your skin. You can apply **Lotion Oxygène** to gain better efficiency. Adjust the vibration intensity based on the size of the area to be cleaned.

LIFT

To revitalize your skin, use the Lift setting. You can use **Sérum La Crème** for additional efficiency.

IPL

Select the IPL mode for glowing, rejuvenated skin. For increased efficiency, use **Sérum La Crème**.

EYE CARE

Select Eye Care to reduce dark circles and wrinkles around the eyes. You can improve your efficiency by using **Sérum La Crème**.

MAINTENANCE

Only use the included accessories to charge the device's battery.


Allow 1 hour for the device to charge.

The indicator light indicates that the battery has been completely charged.


After usage, clean the device's head with the tissue.

CLEAN MODE


- 1** Remove the cotton stopper from the device and cover it with a cotton pad. Replace the cotton stopper to fasten the cotton pad.

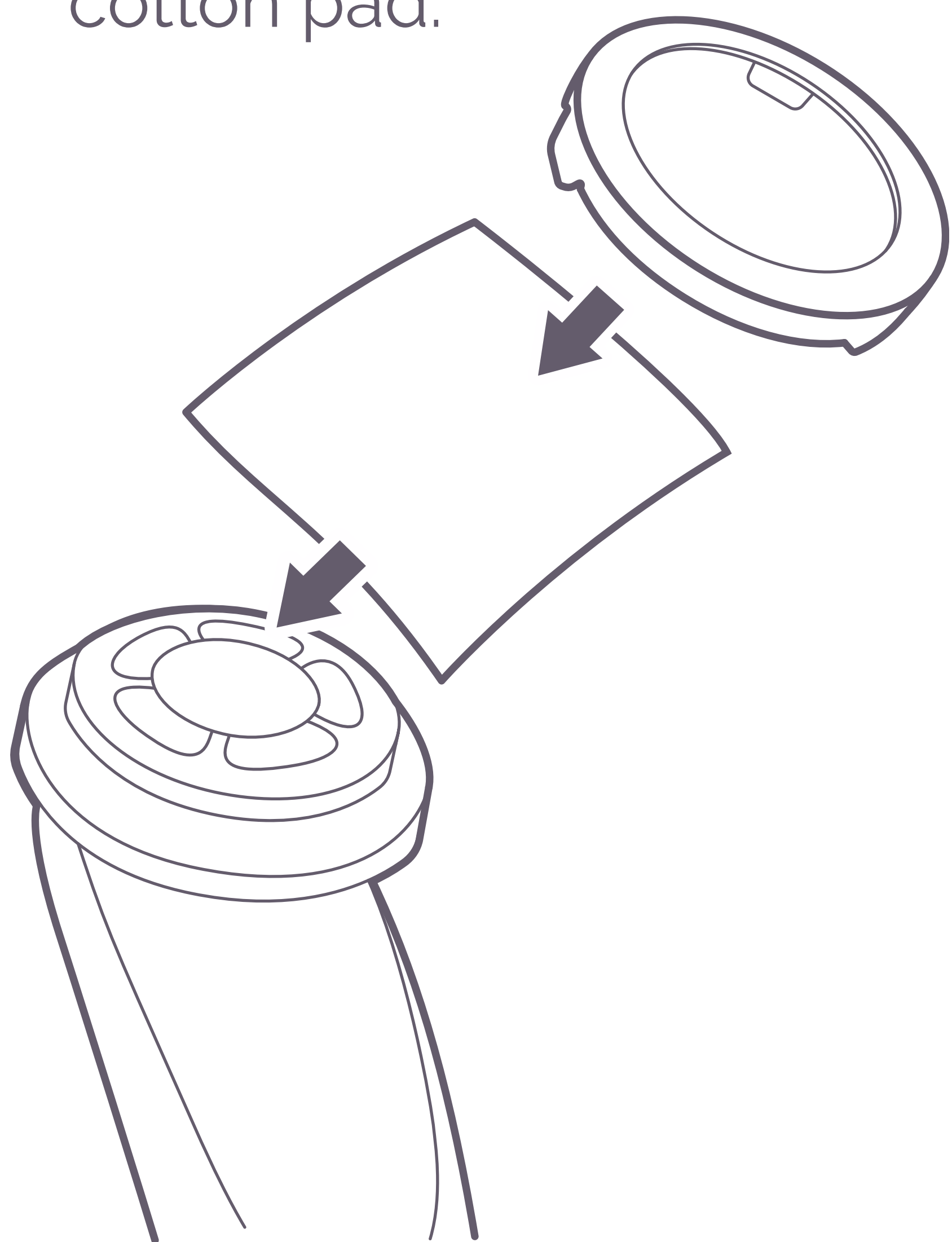
- 2** To improve efficiency, immerse a cotton pad with **Lotion Oxygène**.


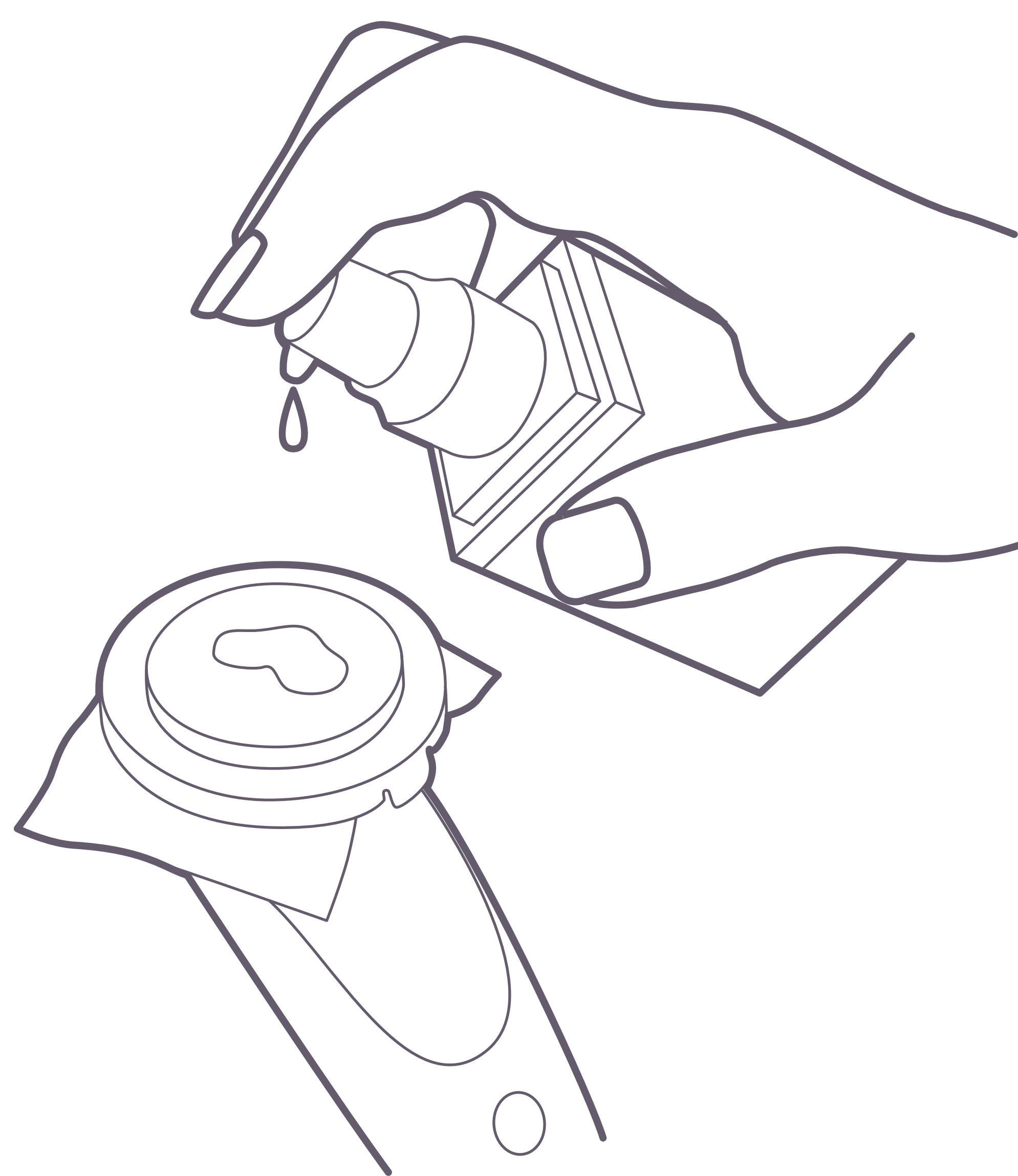
- 3** Choose the **CLEAN** mode by clicking the Power button.

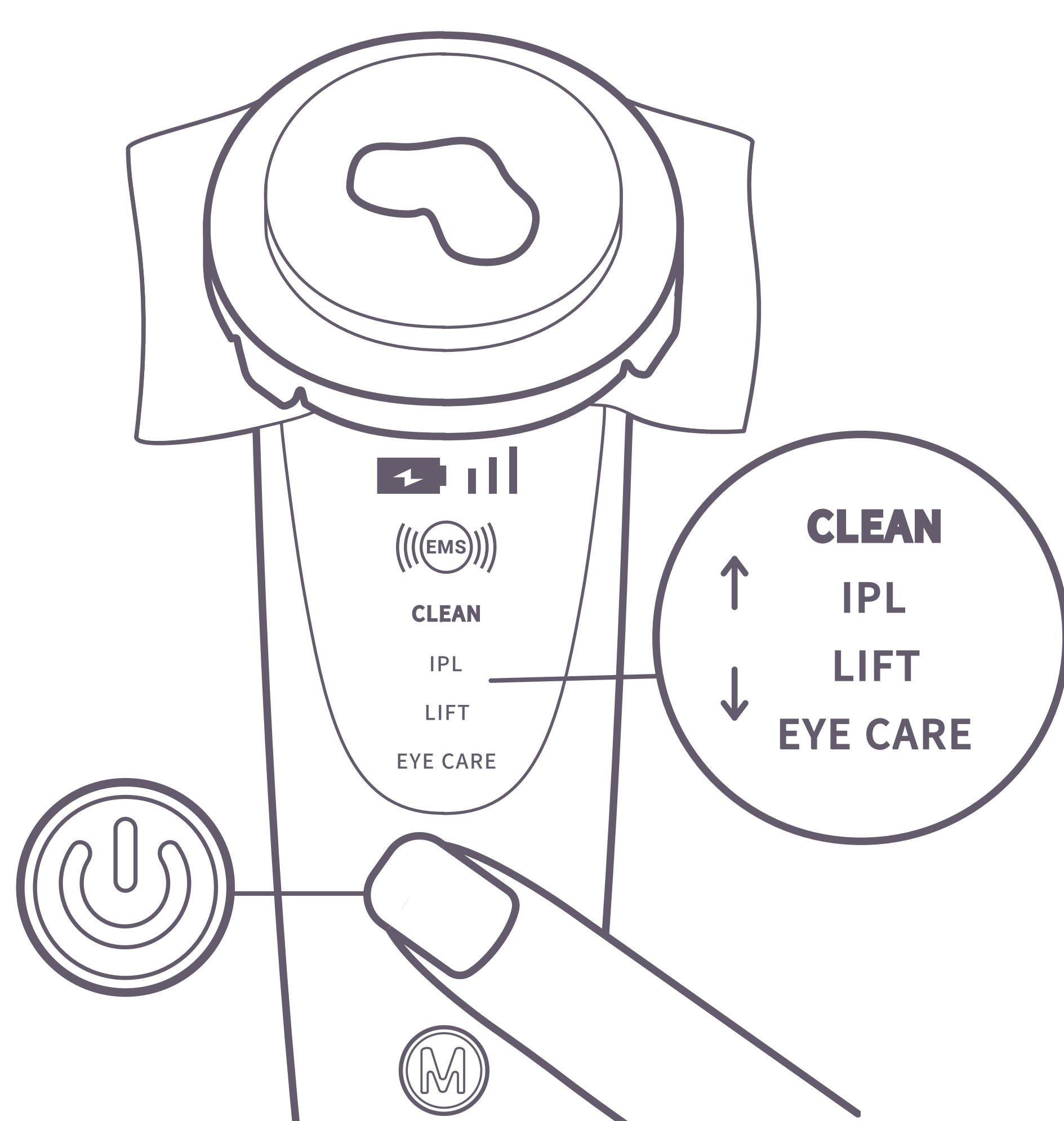- 4** Select vibration level by pressing M.

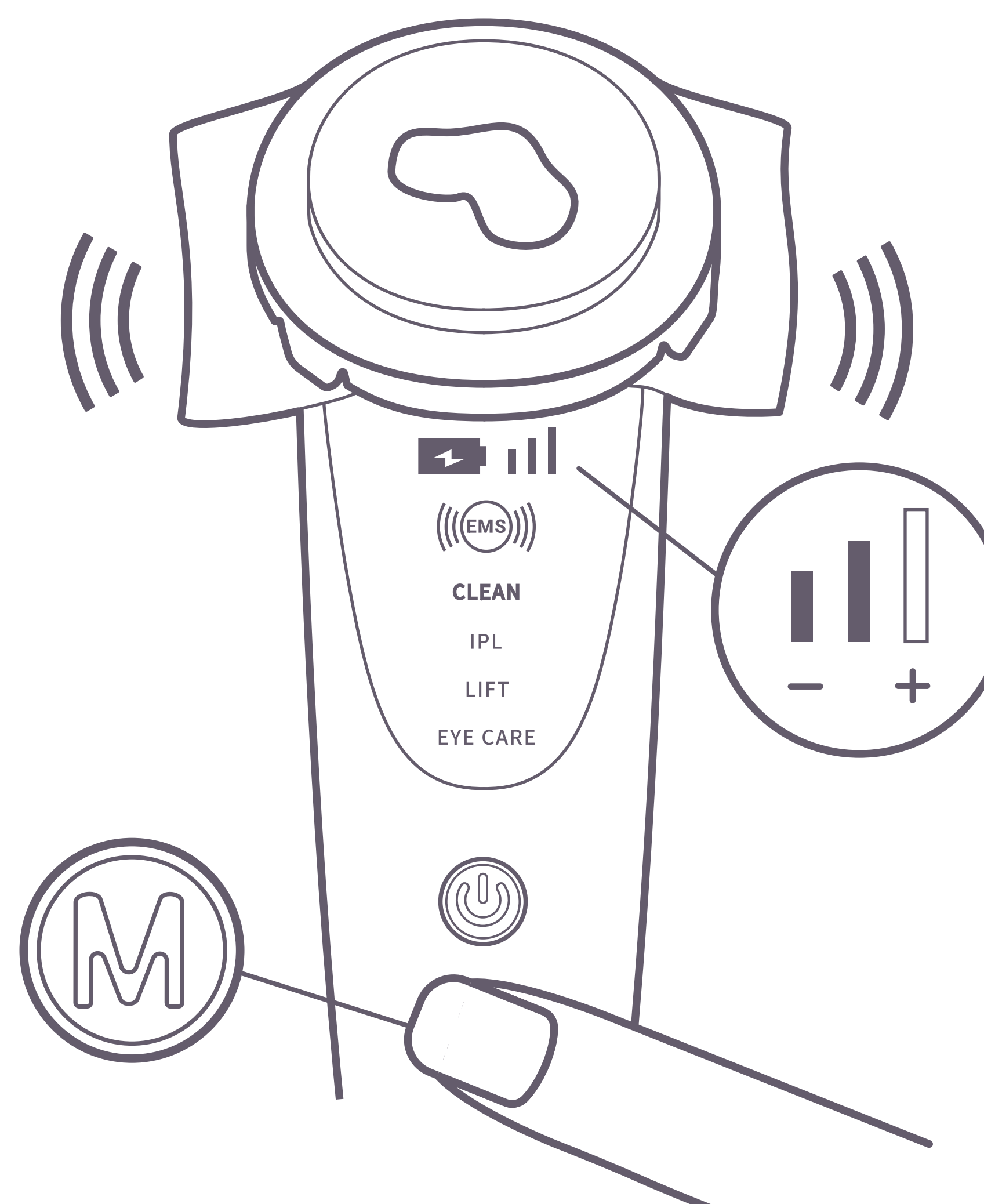
Use Clean mode for one minute to effectively cleanse the skin.

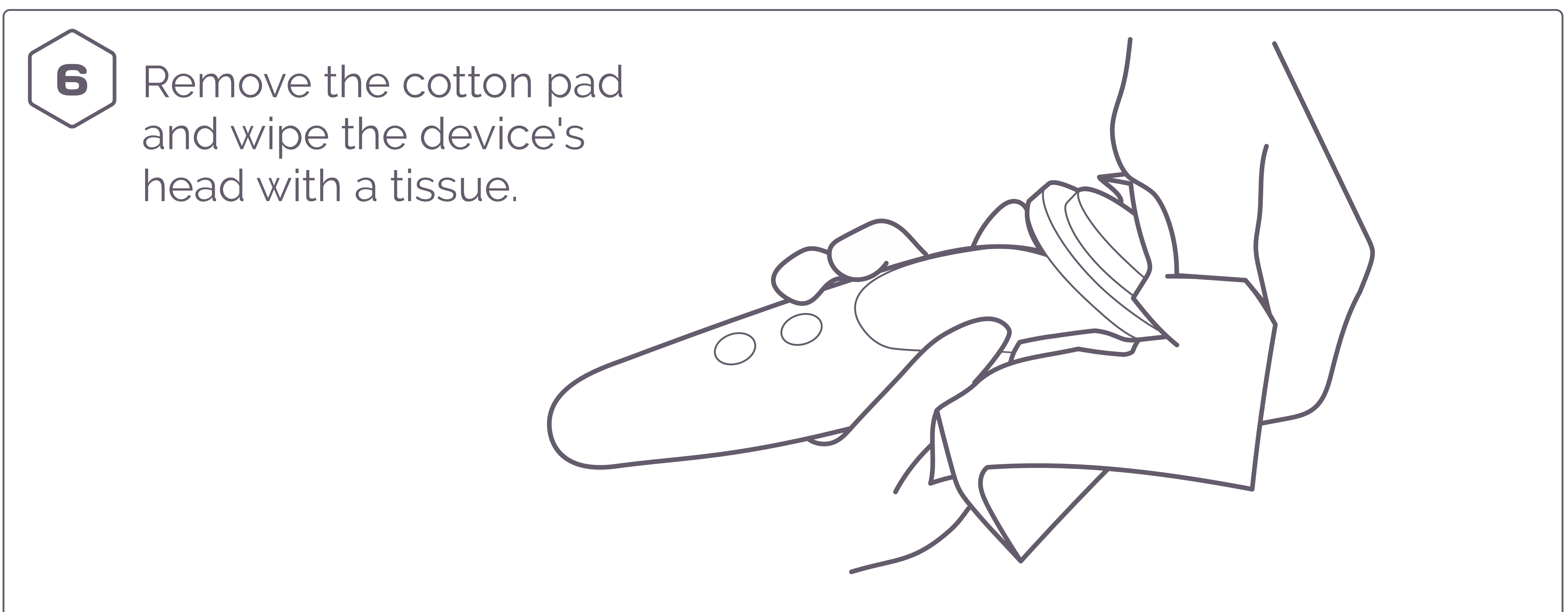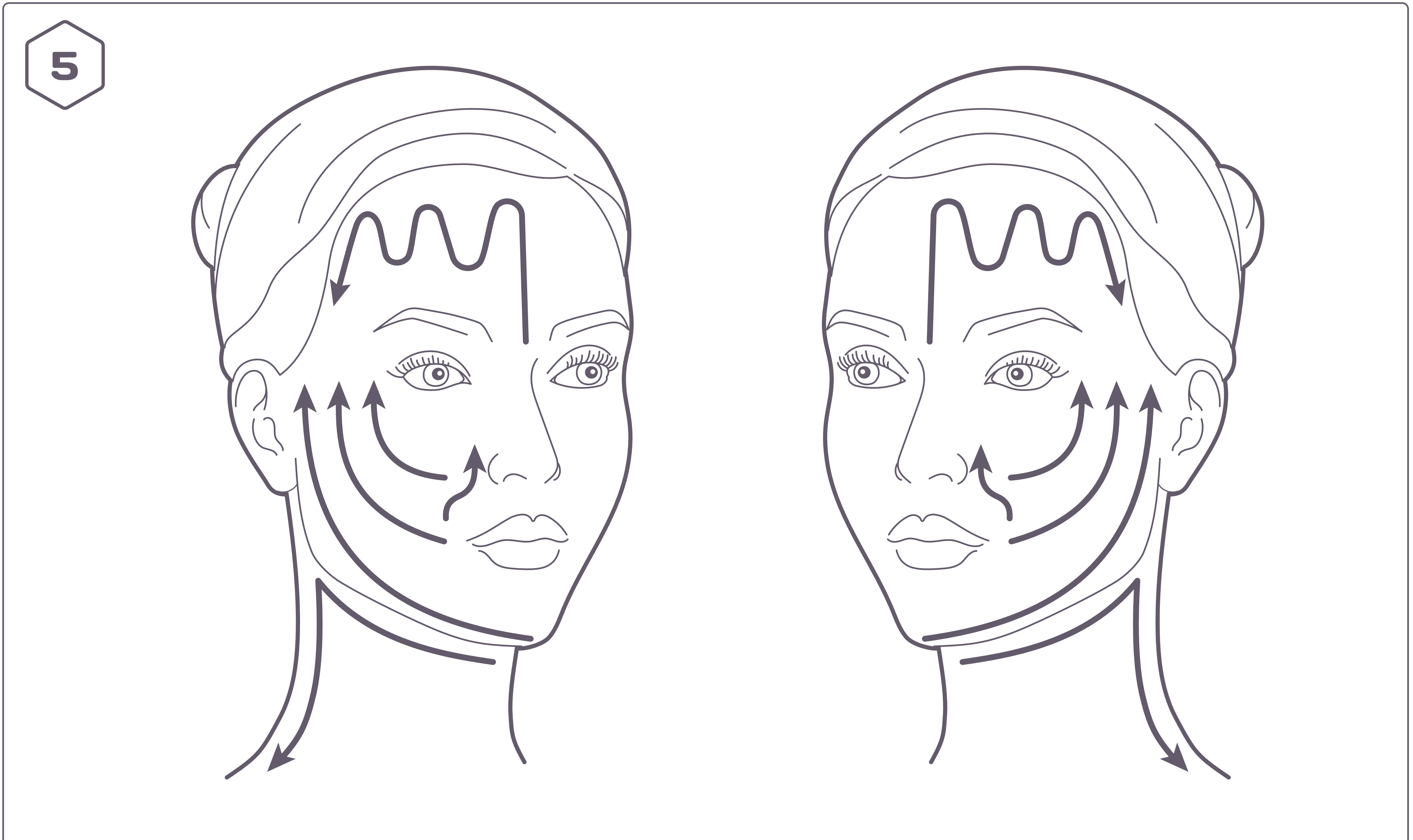
IPL MODE


- 1** Take the cotton stopper from the device and replace it with a fresh cotton pad. Replace the cotton stopper to fasten the cotton pad.

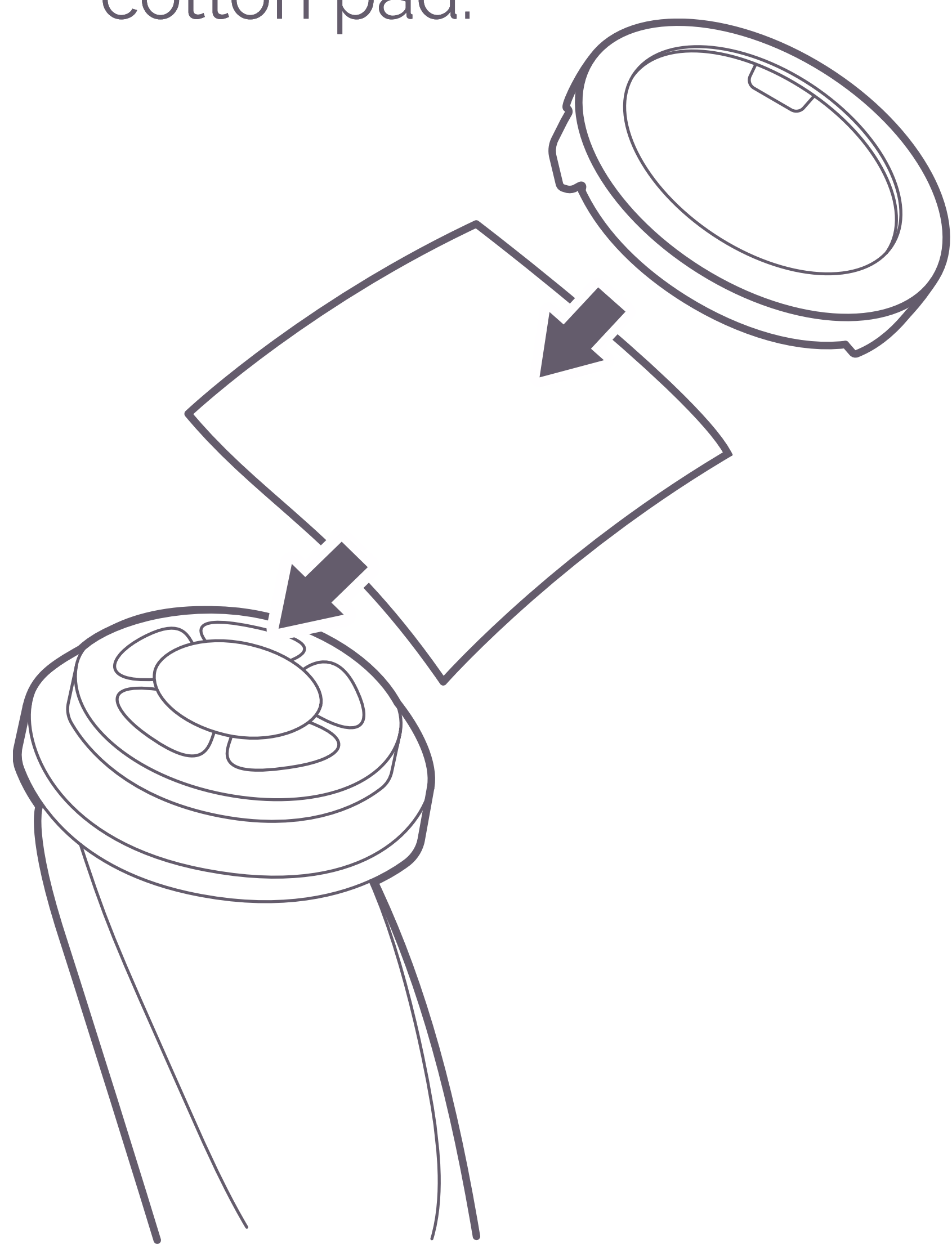
- 2** You can enhance the effectiveness of fresh cotton pads by soaking them with **Sérum La Crème**.

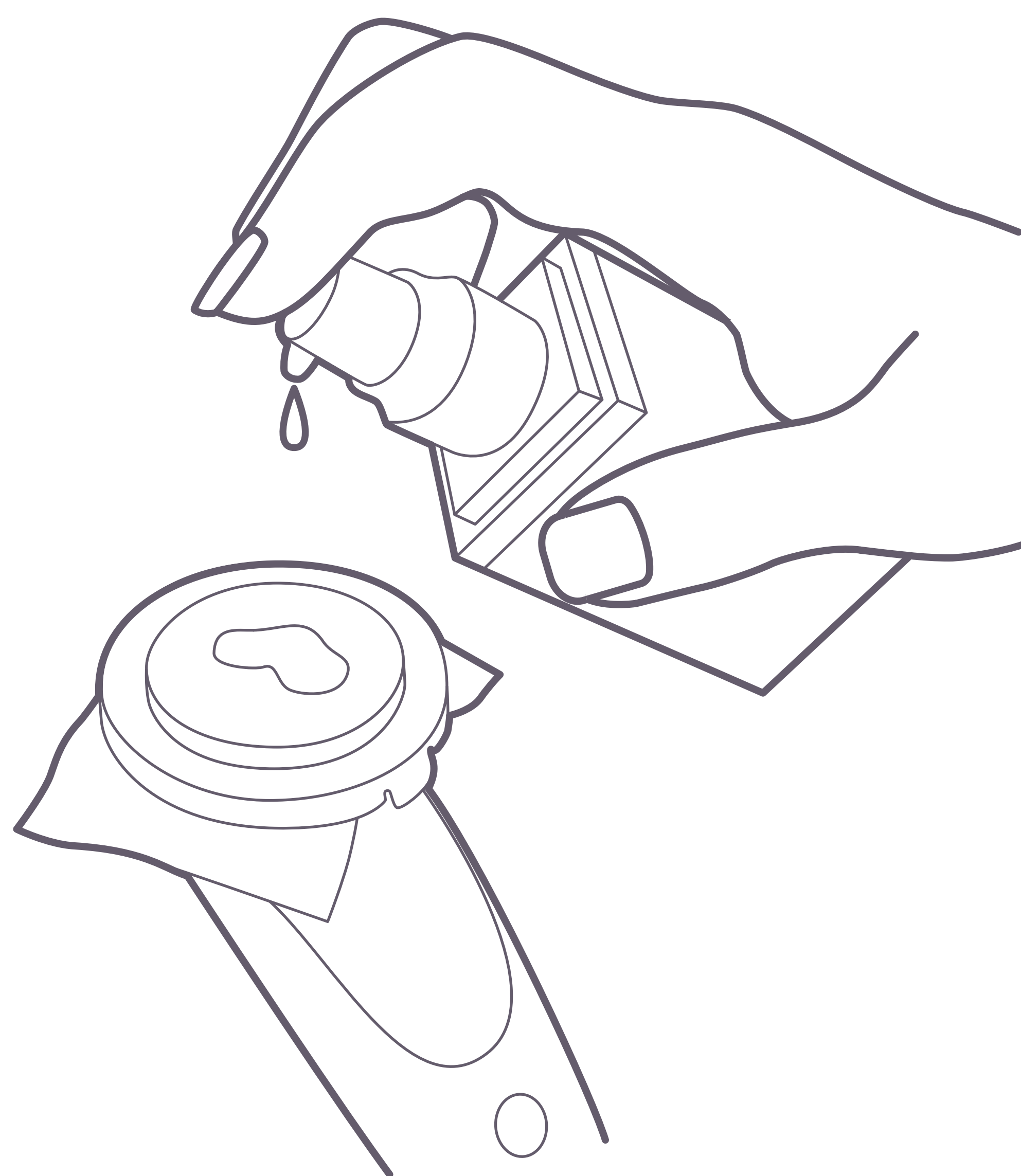- 3** Choose the **IPL** mode by clicking the Power button.

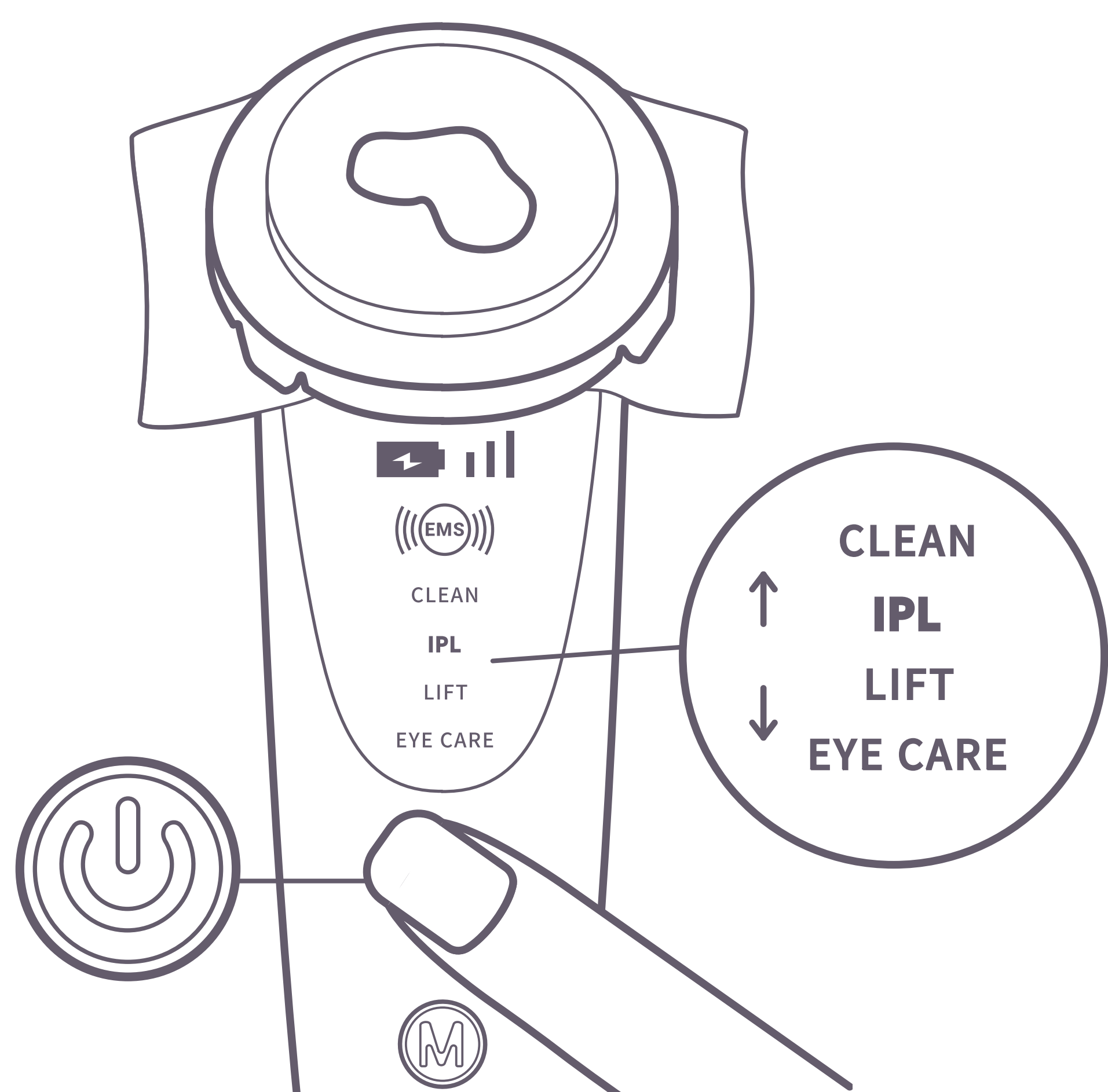
- 4** Select vibration level by pressing M.

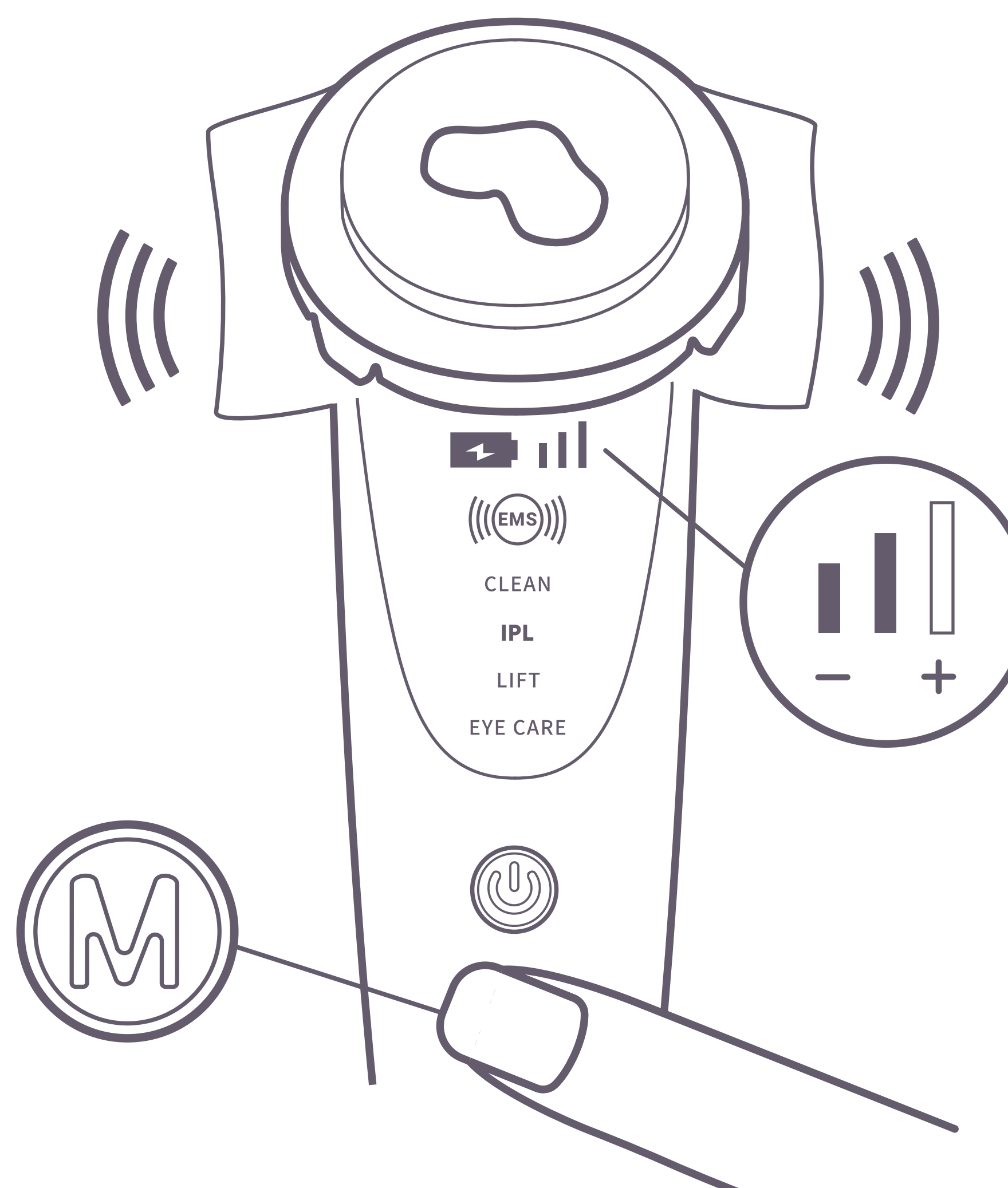For optimal effectiveness, use the IPL mode for two minutes.

- 5** Mode the device sequentially from areas 1 to 6, once for each region.

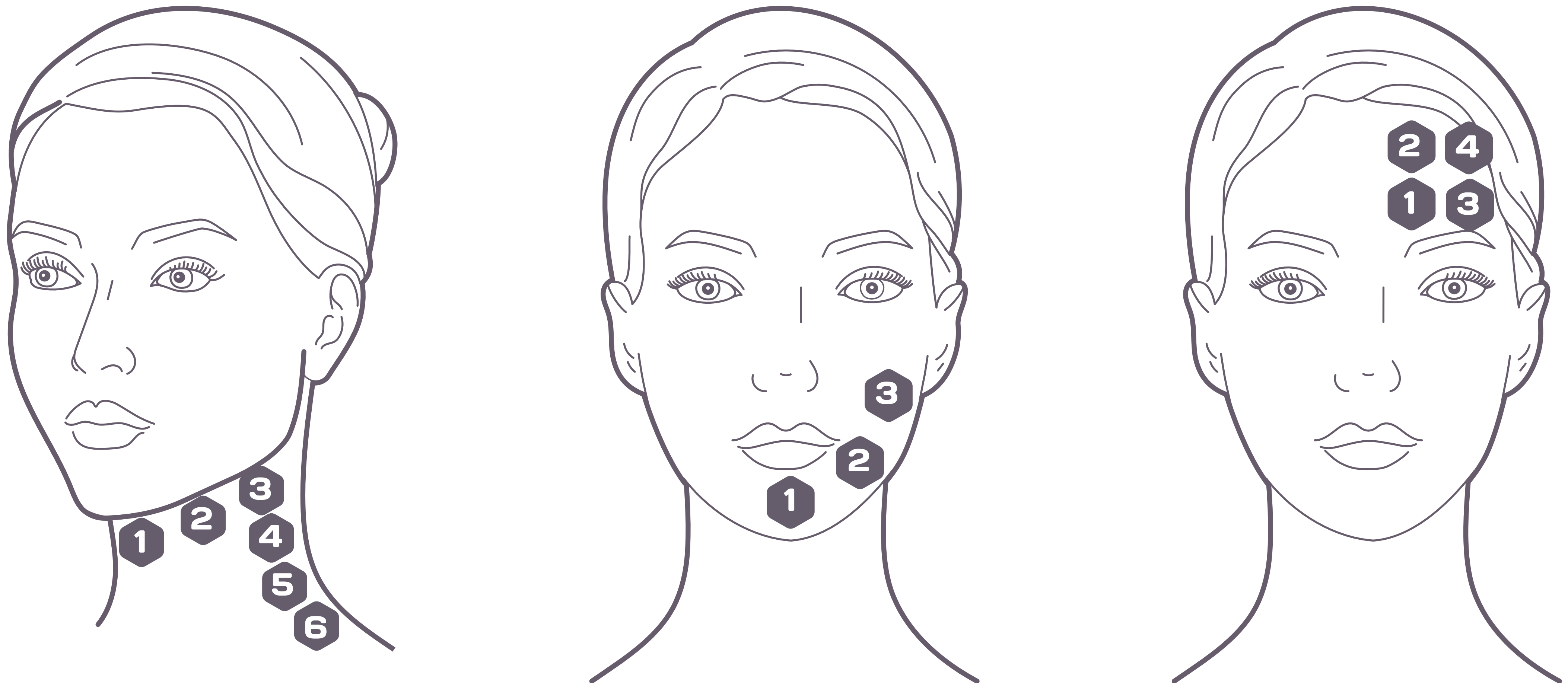
- 6** Remove cotton pad and wipe device's head with tissue paper.


LIFT MODE


1 Remove the cotton stopper from the device and replace it with a new cotton pad.
Replace cotton pad to repair cotton pad.

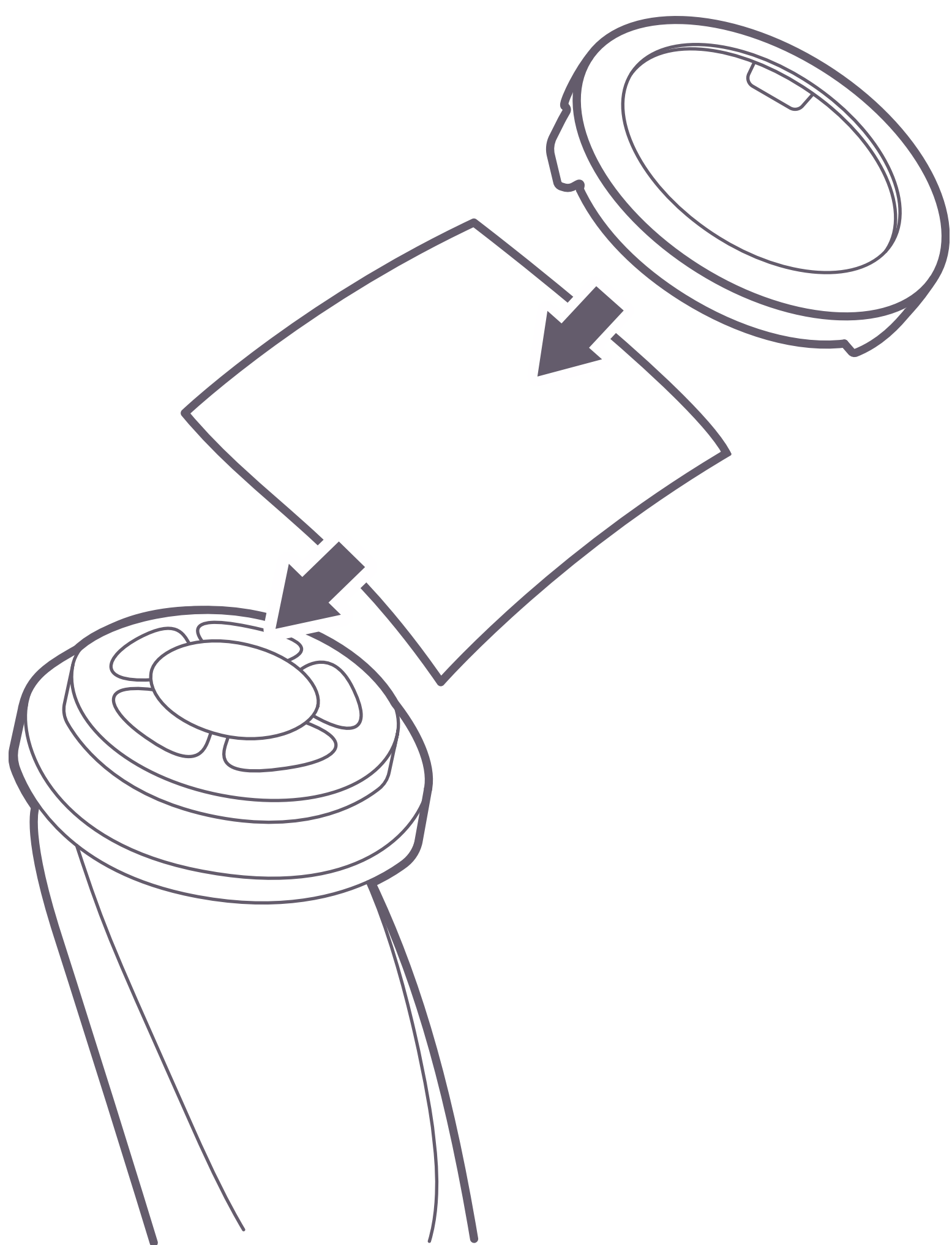
2 New cotton pads can be soaked with **Sérum La Crème** to increase their effectiveness.

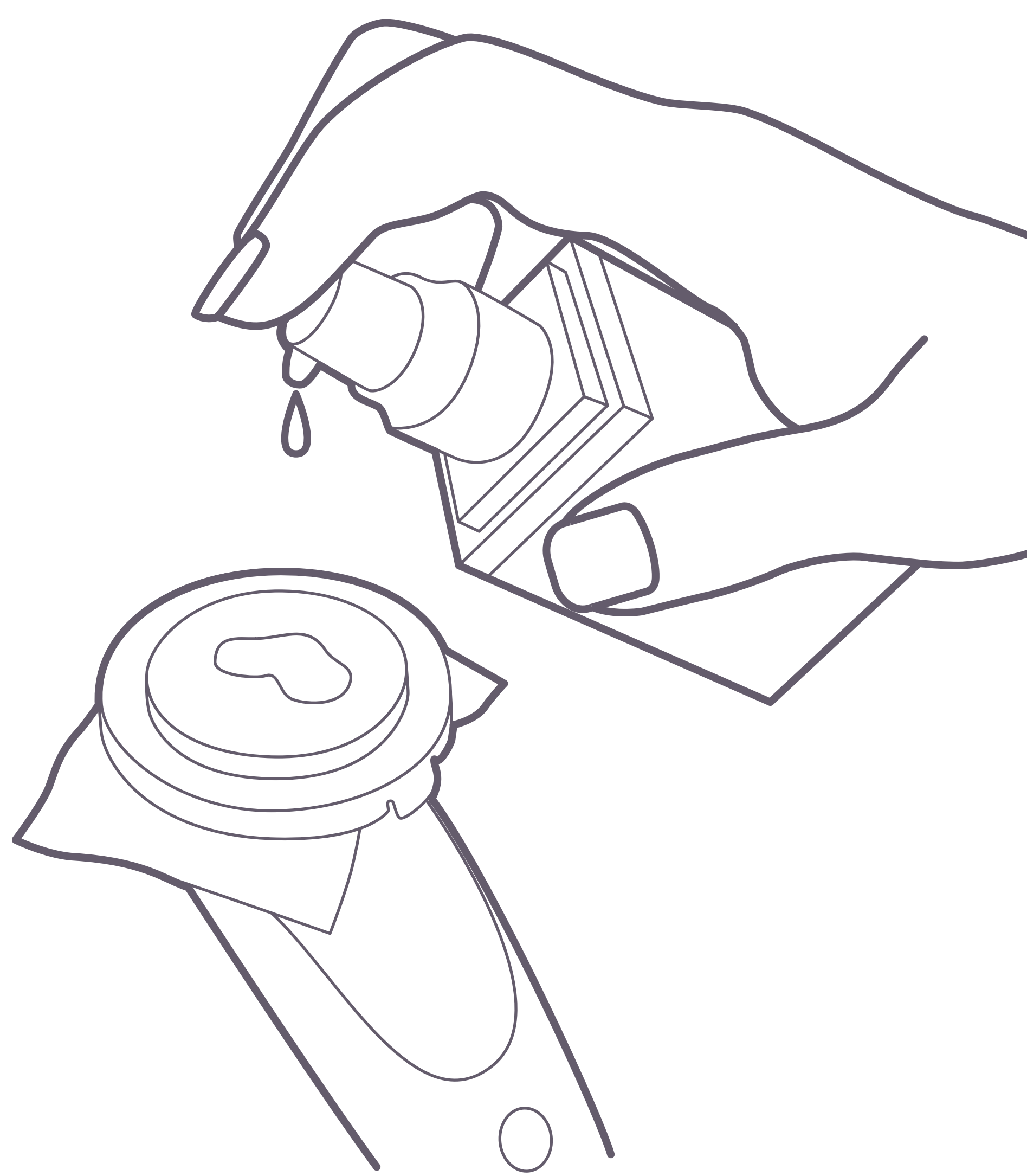3 Select **LIFT** mode by pressing the Power button.

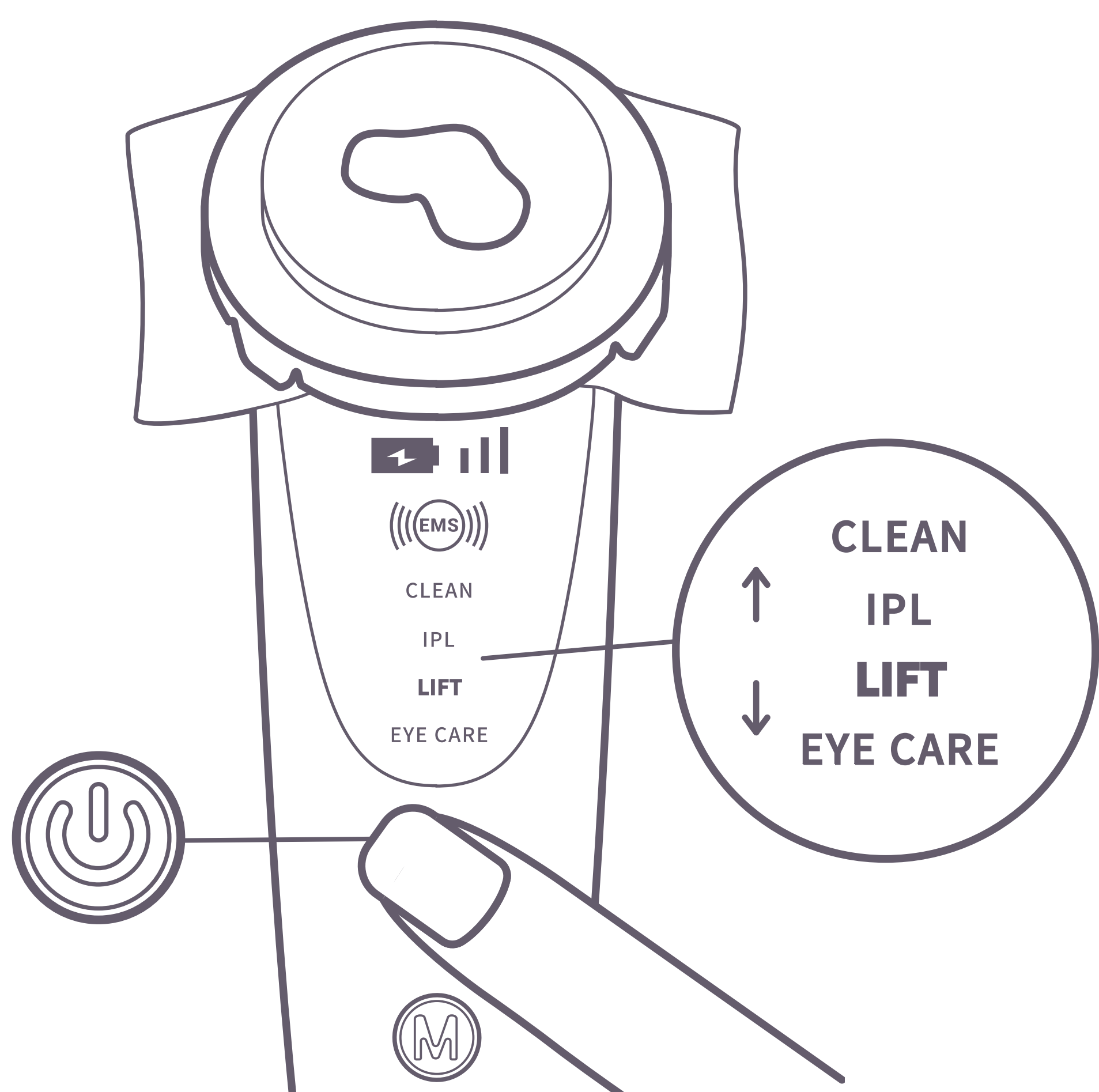
4 Select vibration level by pressing the M button.

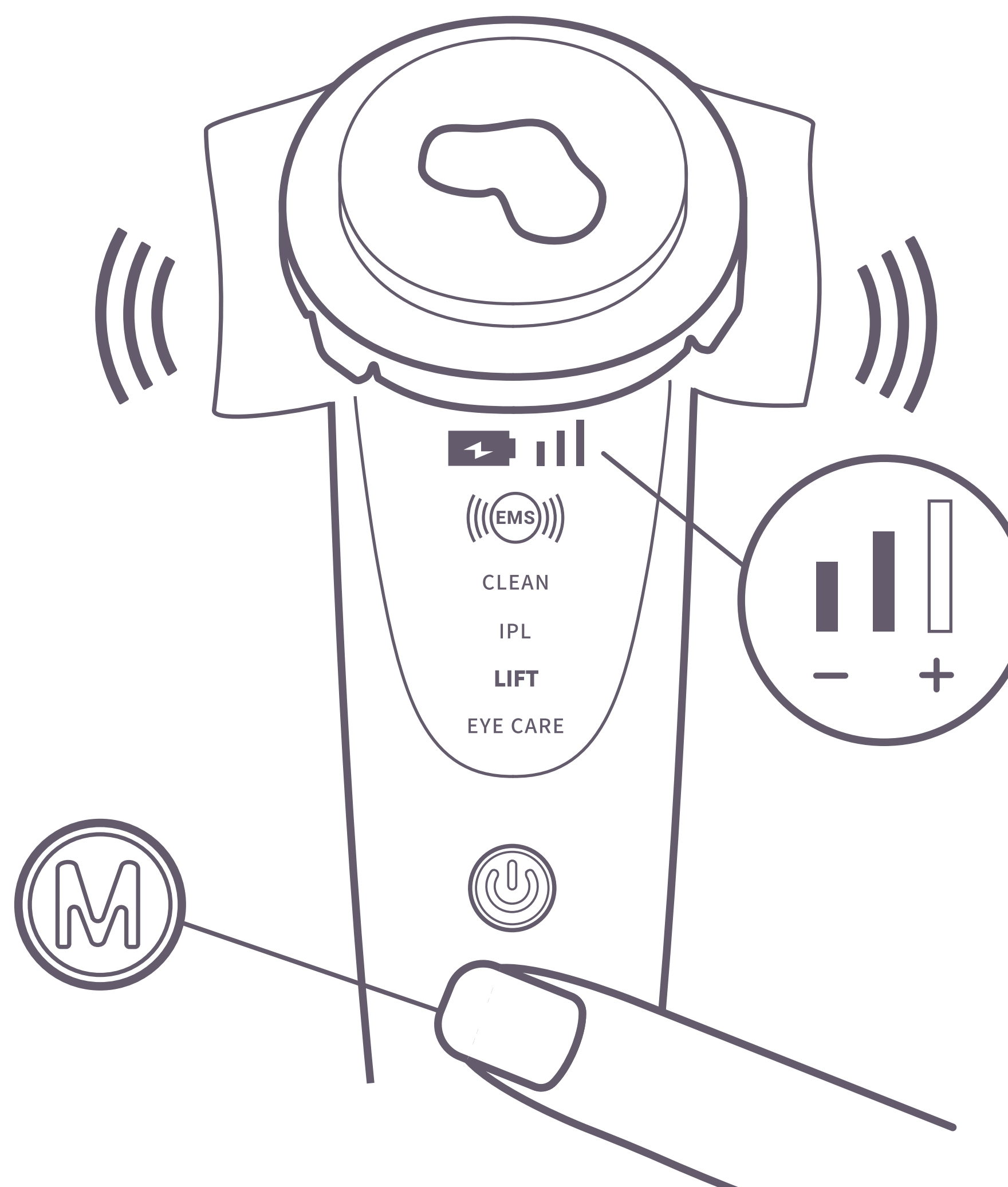Use Lift mode for 3 minutes for efficient lifting.

5 Mode the device sequentially.

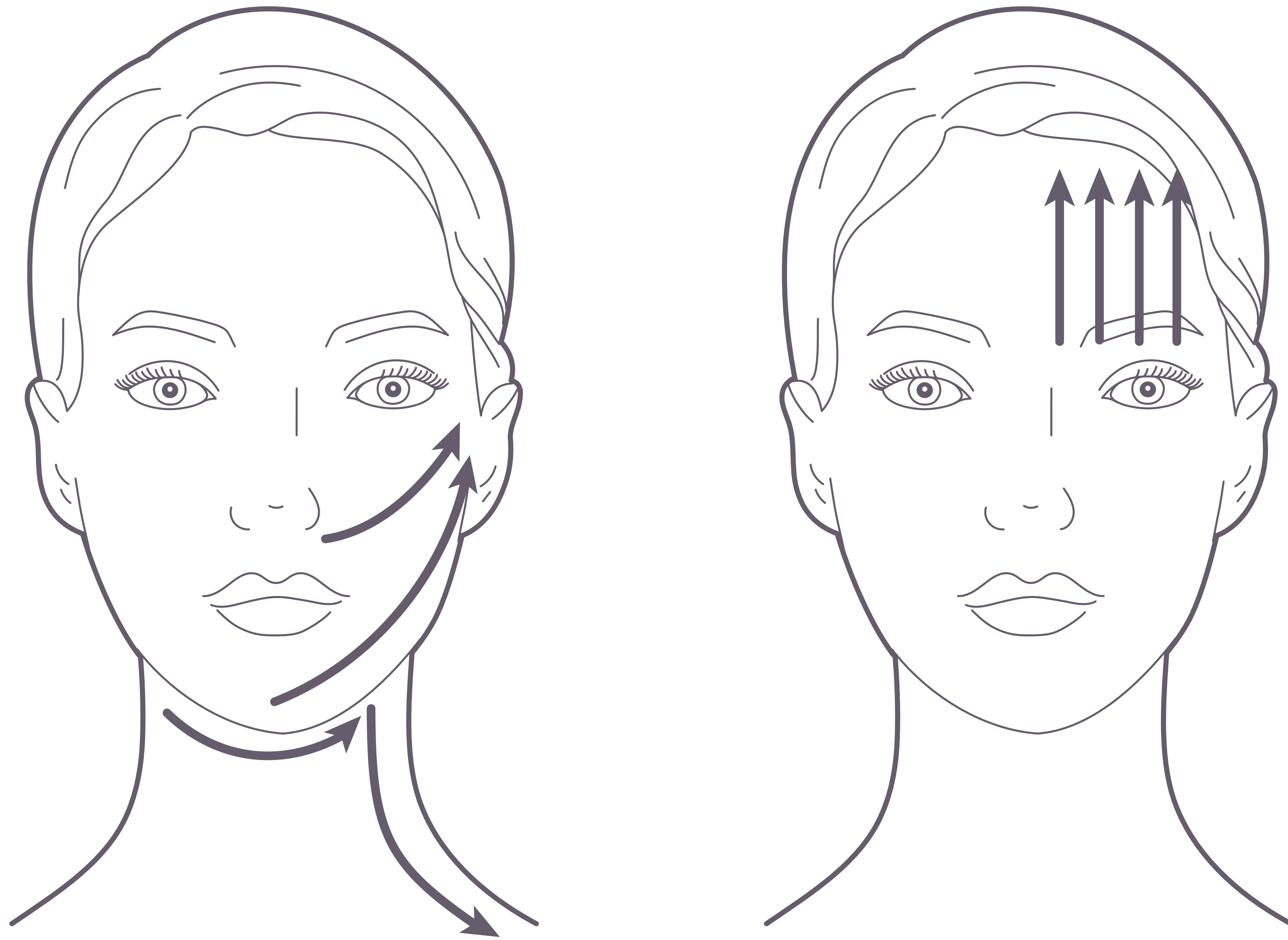
6 Clean the device's head with a tissue after removing the cotton pad.


EYE CARE MODE


- 1** Take the cotton stopper from the device and replace it with a fresh cotton pad. Replace the cotton stopper to fasten the cotton pad.

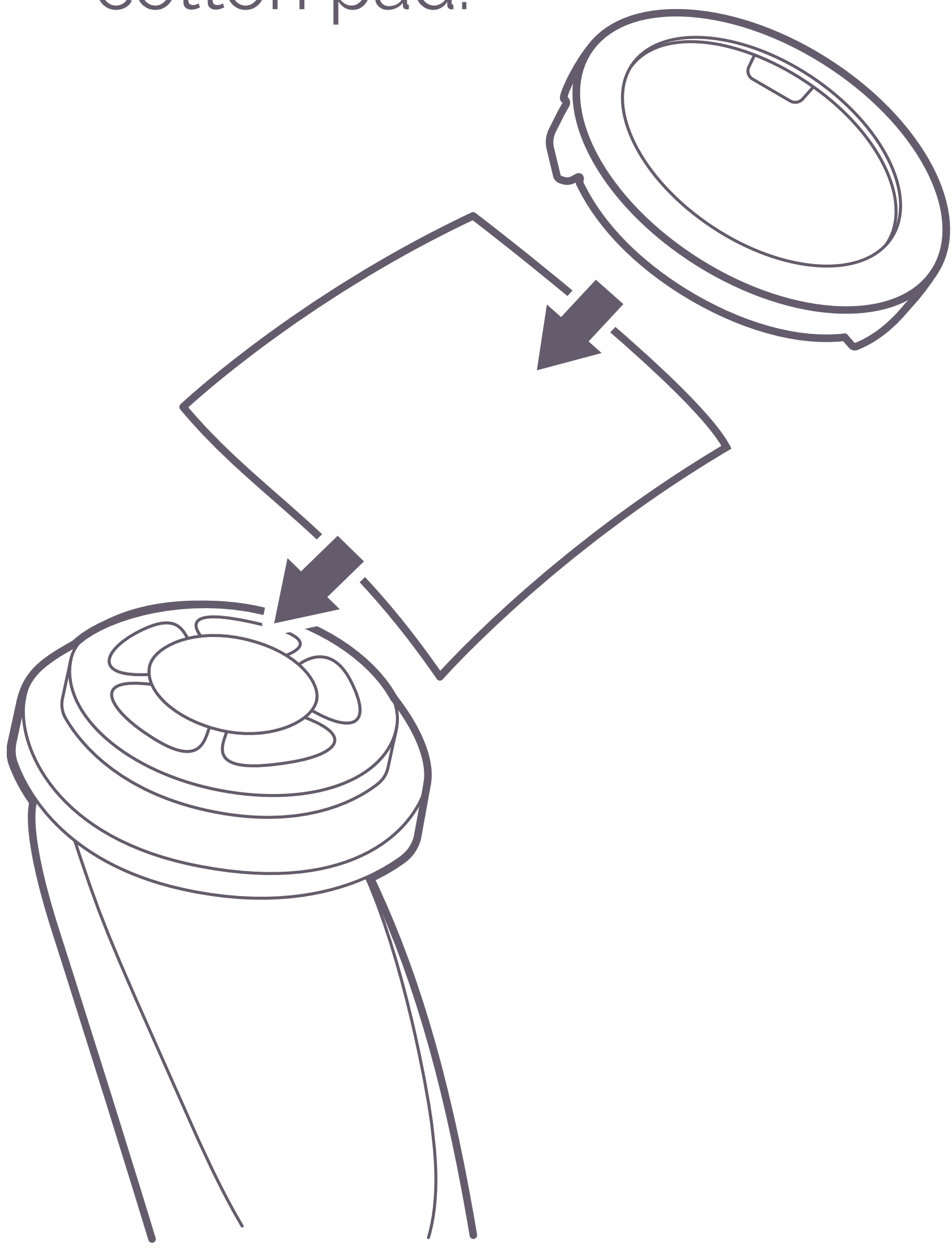
- 2** You can enhance the effectiveness of fresh cotton pads by soaking them with **Sérum La Crème**.

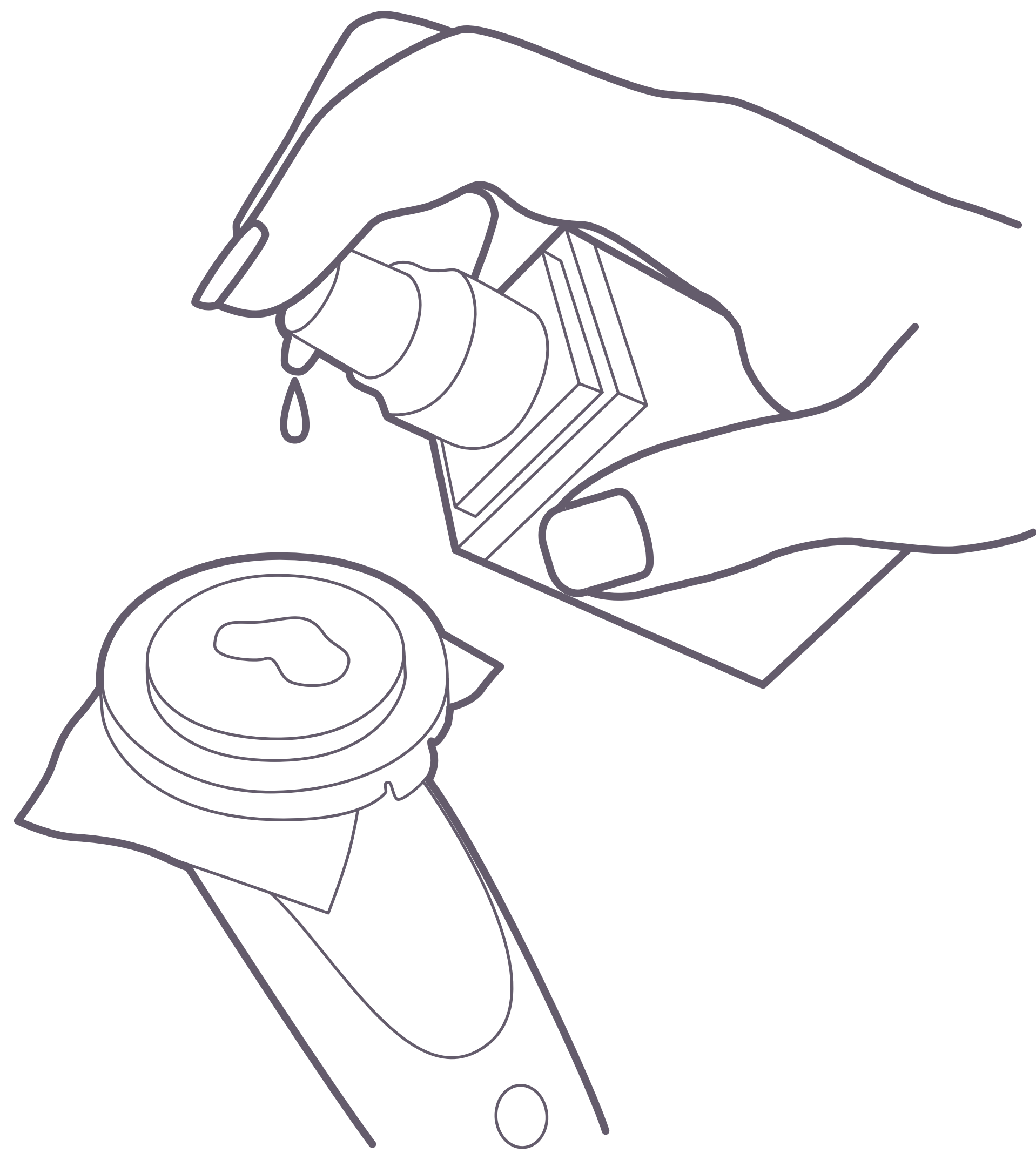- 3** Select **EYE CARE** mode using the Power button.

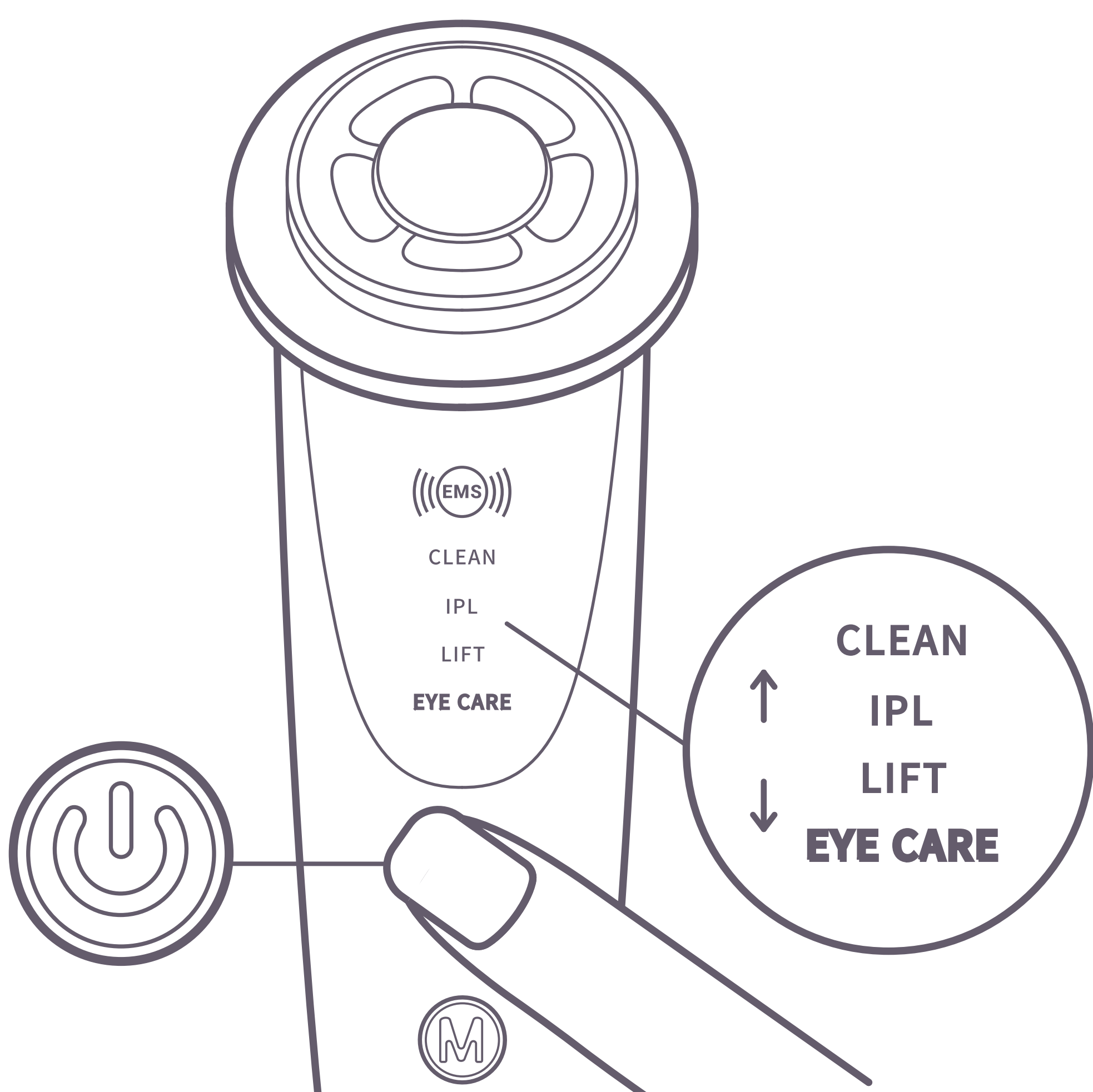
- 4** Select whether or not to engage the red-light therapy by pushing the "M" button.

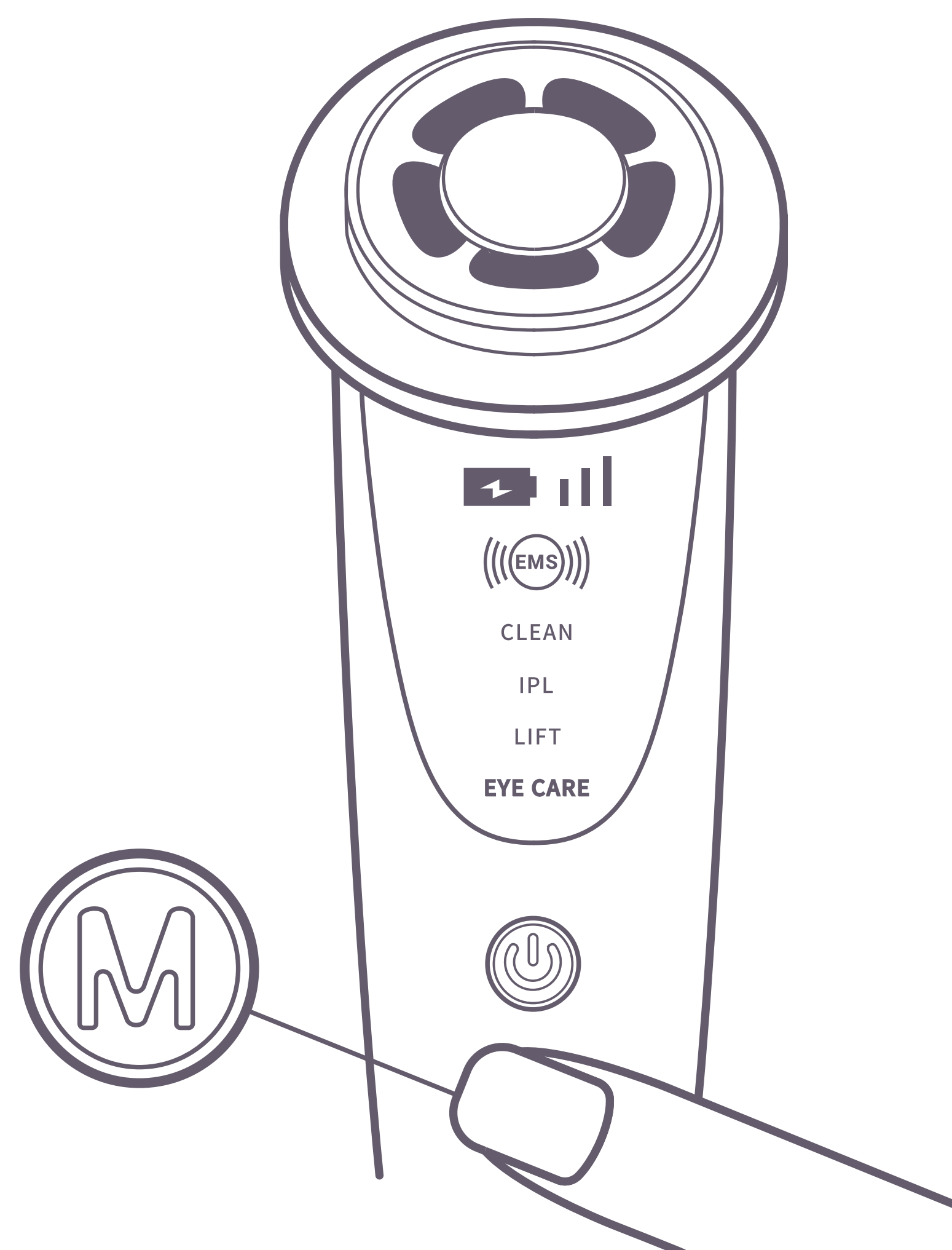
Use the Eye care mode for 2 minutes for maximum efficiency.

5 Treat each region at the same time.

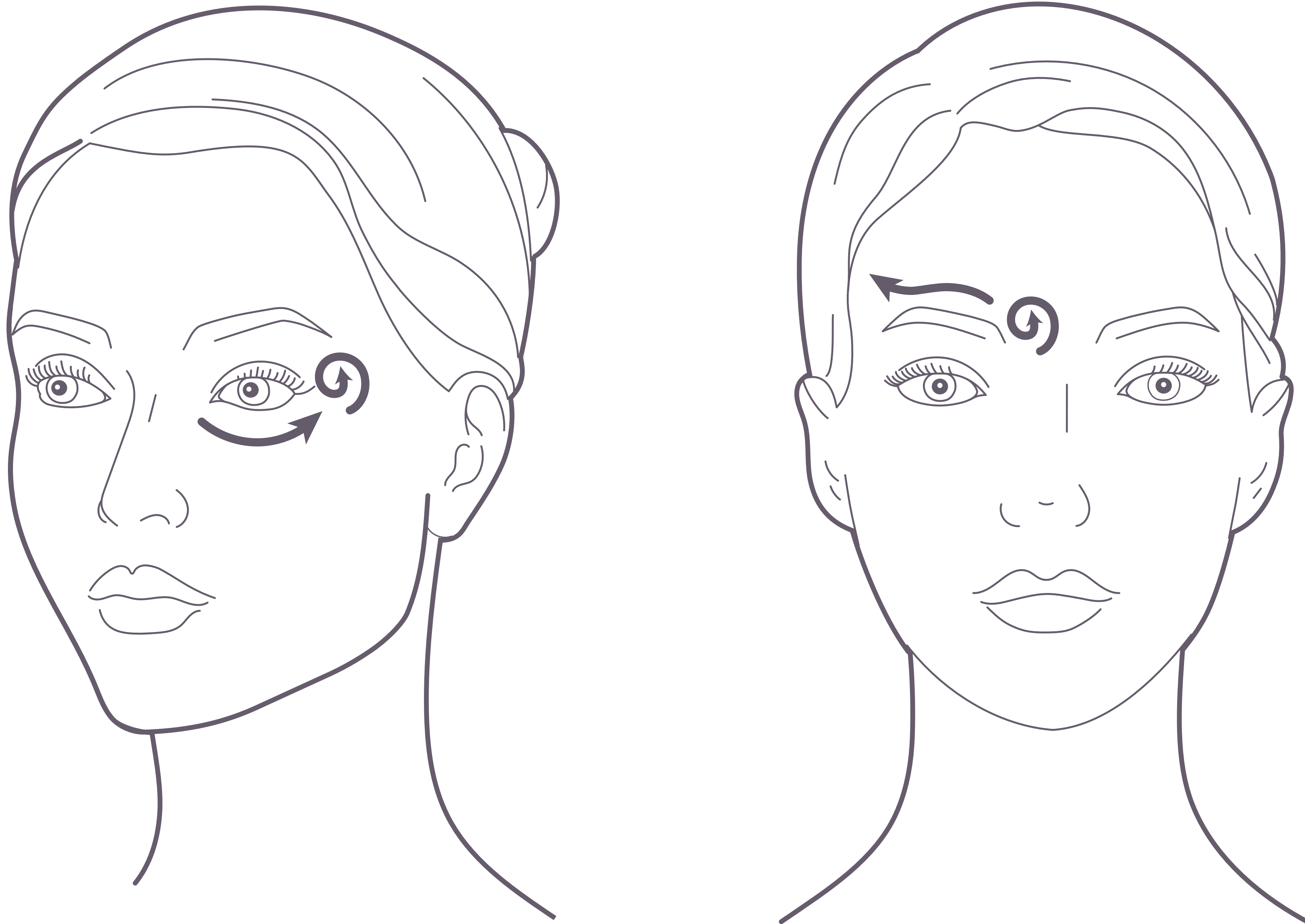6 If you used a cotton pad, remove it and wipe the device's head with tissue.

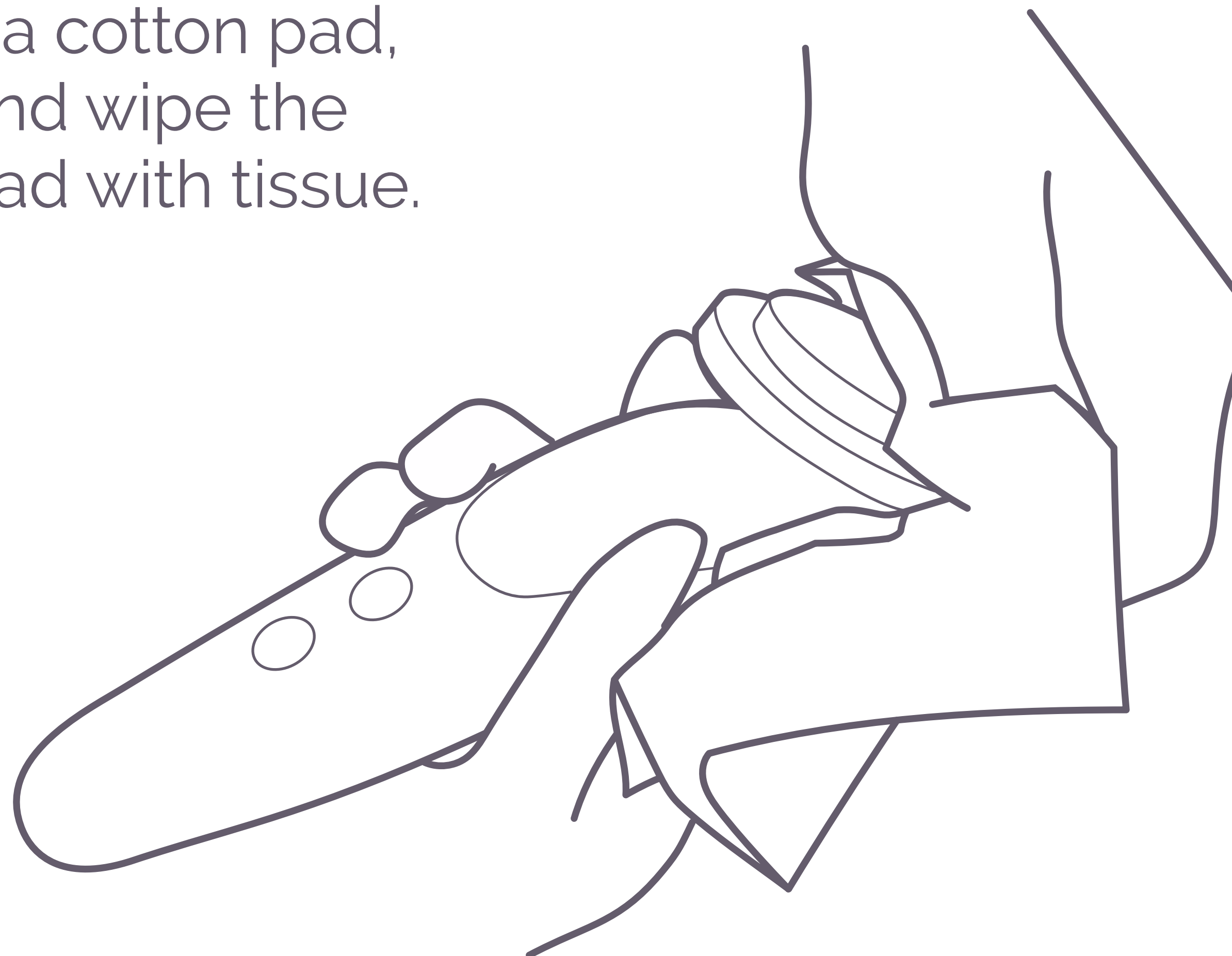SAFETY PRECAUTIONS

DO NOT USE THIS UNDER THE FOLLOWING CONDITIONS:

- In conjunction with pacemakers or other implanted electronic medical devices, artificial heart-lung machines, life-support medical devices, electrocardiographs, or other wearable gadgets.
- If you have a skin ailment or problem, an allergy, or sensitive skin.
- If you're ill, inebriated, pregnant, menstruating, breastfeeding, or taking a drowsy medication.
- If you've had eye or face surgery or are being treated for eye or face diseases.

DO NOT REPEATEDLY APPLY TO THE SAME AREA:

- Ensure that the device is being moved across the skin while in use. Do not use on areas that cannot be reached by the entire surface of the electrode head.

WARNING!

- Do not use near combustible materials.
- Do not use in the bathroom or near a bathtub, shower, sink, or other water-filled containers, as it is not water-resistant.
- Do not attempt to disassemble, repair, or modify this device.
- Stop using if any problem occurs, such as an abnormal odor or excessive heat.
- If you feel an earthquake or witness lightning while a device is charging, switch off device.
- Store the equipment in a suitable location.

ATTENTION!

- Ensure the electrode head is not cracked, distorted, or fractured before usage.
- Continuous use should not exceed eight minutes.
- Prevent the device from coming in touch with metal surfaces, other beauty products, or electronic appliances while the power is on.
- Remove your makeup before using the device on your face.
- At least 24 hours should elapse between subsequent applications to the same region.
- Before cleaning the device, ensure the power is off, and the AC adapter is unplugged.

POWER SUPPLY

- Use the provided AC adapter and USB cord with this product.
- When storing the device, do not wrap the USB cord around the device or AC adapter.
- Before charging the device, ensure that the charging port and USB cord are dry.
- Do not pull the AC adapter or USB cable when moving or lifting the device.
- Ensure operation within the voltage range of 100 to 240 volts.

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